BÁSCULA INALÁMBRICA DE ANÁLISIS CORPORAL  
WIRELESS BODY FAT SCALE  
BALANÇA DE ANÁLISES CORPORAL SEM FIO  
ÉCHELLE D'ANALYSE CORPORELLE SANS FIL

## IVSCB

- ES** GUÍA DE USUARIO ..... 2
- EN** USER GUIDE ..... 7
- PT** MANUAL DE INSTRUÇÕES..... 11
- FR** GUIDE DE L'UTILISATEUR..... 15

## LA CAJA CONTIENE:

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## SOME TECHNICAL DETAILS

**Dimensions:** 305 x 305 x 24 mm

**Bluetooth Smart 4.0 communication protocol with the IvyHealth application**

**Fitted with high precision sensors**

**Limitless registration of users in the IvyHealth application**

**Weight measurement range: 3-150 Kg**

**Body fat measurement range: 5% - 45%**

**Body water percentage range: 29%-69.5%**

**Bone mass/muscle mass index: 0.1Kg**

**LED screen and "low battery" indicator**

**It works on 4 AAA batteries (not included)**

**Automatic reset to "0"**

**It displays the body weight reading on the screen and readings of the other parameters on your mobile phone or tablet with the IvyHealth application.**

## INTENDED USE

IVY HEALTH's IVSCB body analysis scale was designed for use by adults. With its non-invasive electronic technology, it offers a series of body values which are displayed on the scale itself - weight - and others on the free IvyHealth application. (See further on)

## SAFETY WARNINGS

1. Your feet must be dry and bare when using the scale for more precise measurement.
2. This scale has high precision sensors. The scale must be placed on a smooth, rigid and level surface for it to work properly. For example, a tiled floor is ideal, whereas a deep pile carpet or rug could distort some of the scale's measurements.
3. Do not move the scale while it is switched on. It could cause measurement errors.
4. Do not allow children or the infirm, who could fall off it, to use the scale without your help.

## VERY IMPORTANT


1. You are not recommended to use the scale if you are pregnant.
2. This scale emits an imperceptible and harmless electrical current in order to make the advanced measurements. Nevertheless, its use is not recommended if you use a pacemaker or have any electronic implant in your body. It could affect its proper working.

3. Don't move the scale while it is activated. It could cause measurement mistakes.
4. Don't let the children or people who can fall off the scale use the device without your supervision.
5. If *LO* appears on the screen, change the batteries. There are wearing out and the results may not be so accurate.

### Cleaning the device:

1. Use a damp cloth with a soap solution, but make sure no liquid gets into the device.
2. Do not use abrasive liquids or cloths or hard utensils when cleaning.
3. Do not expose it directly to the sun, extreme humidity and/or dust.

### LED DIAGRAM

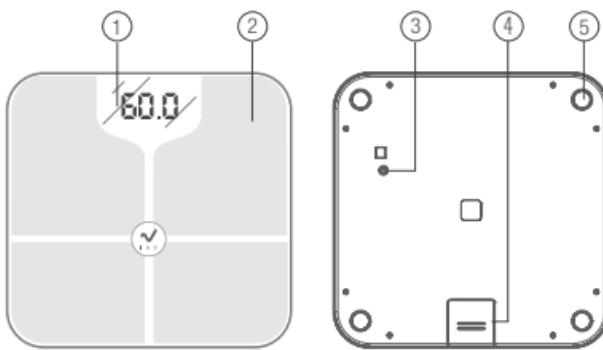- ① LED Screen  
② Conductive areas  
③ Change measurement unit  
(Kilograms, Pounds, Stones)  
④ Battery compartment  
⑤ Stabilizers

### WHAT CAN I DO WITH THIS SCALE?

For starters, you can weigh yourself with this scale. Open the lid of the battery compartment and place 4 AAA batteries in the position indicated. Simple: the scale is ready to use, just like any normal scale.

However, the reason you bought this scale is that you are looking for something more.

You can use the scale without the IvyHealth app, but it will only tell you your weight and body mass index.

Download the free **IvyHealth** application from the Apple or Google Play Store. When you first open the application, you must register and set up your user name.

It is very easy if you follow the instructions on your mobile phone's screen.

The IvyHealth application can be automatically paired with your recently acquired smart scale.

1. To make a measurement, press the icon on the main screen. You will access the measurement screen and can now step on the scale with your bare feet.
2. The scale screen displays your weight. While your weight flashes on the screen, the scale continues to gather information on the other values it can offer.

3. Once the transmission has been completed, you can see the following data on your mobile phone screen:
  - a. Your weight
  - b. Your BMI
  - c. Your muscle mass
  - d. Your body water percentage
  - e. Your visceral fat
  - f. Your body fat
  - g. Your bone mass
4. Press to save and your first reading is recorded for the future. You can consult your phone whenever you wish. You will be protected by your own password.
5. Subsequent measurements will be kept in your mobile phone whenever you weigh yourself with the application opened and you follow the previous steps.

### **What else can I do with the application?**

If you use the scale together with the free **IvyHealth** application, you can obtain the readings for all the body values listed in the previous paragraph. However, the application also allows you to:

- 1. Consult your statistics.** Just click on the historic measurements section.
- 2. Save and display the results of as many users as you wish.** Every time somebody makes a measurement from their mobile phone, and if they are registered as a user, they can record their results and protect them with a password

## In addition to the scale

The IVY HEALTH package of health complements is made up of different devices affording you greater control over your health.

Check out our latest products for  
healt care in this link:  
<http://www.ivyhealth.es>

## LEGAL NOTICE

### **Declaration of Conformity**

FONEXION SPAIN S.A. declares that the IVSCB wireless smart scale is in compliance with the essential requirements and other relevant provisions of Directive R&TTE 1999/5/EC, Directive EMC 2004/108/CEE and the directive on the restriction of the use of certain hazardous substances in electrical equipment, (ROHS) 2011/65/EU. You can download a copy of the complete declaration at: <http://www.fonexion.com/doc>.

### **Guarantee**

The IVSCB wireless smart scale has a two-year guarantee from the date on the sales receipt. The sales receipt is your guarantee, so keep it. Nevertheless, this guarantee does not cover any fault caused by accidents, improper use, wear, negligence, lightning, modifications made to the equipment, or any attempt at adjustment or repair in the home or by unauthorised agents.

### **Bluetooth**

The Bluetooth brand and its logos are registered by Bluetooth SIG, Inc. and use of said brands by IVY HEALTH is made under licence. Other brands and trade names belong to their respective owners.

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Distribuido por **FONEXION SPAIN S.A.** VAT: ESA48937239